

ATS5004D

User manual



ATTENTION!

Measuring directly on the **line voltage** can be very dangerous.

Copyright ©2015 TiePie engineering.
All rights reserved.

Revision 2.13, September 2015

Despite the care taken for the compilation of this user manual, TiePie engineering can not be held responsible for any damage resulting from errors that may appear in this manual.

Contents

1	Safety	1
2	Declaration of conformity	3
3	Introduction	5
3.1	Differential input	5
3.1.1	Differential attenuators	7
3.1.2	Differential test lead	9
3.2	Sampling	9
3.3	Sample frequency	10
3.3.1	Aliasing	11
3.4	Digitizing	13
3.5	Signal coupling	13
4	Driver installation	15
4.1	Introduction	15
4.2	Where to find the driver setup	15
4.3	Executing the installation utility	15
5	Hardware installation	21
5.1	Power the instrument	21
5.1.1	External power	21
5.2	Connect the instrument to the computer	22
5.2.1	Found New Hardware Wizard	22
5.3	Plug into a different USB port	25
6	Front panel	27
6.1	Channel input connectors	27
6.2	Power indicator	27
7	Rear panel	29
7.1	Power	29
7.1.1	USB power cable	30
7.1.2	Power adapter	31
7.2	USB	31
7.3	Extension Connector	31
8	Specifications	33

8.1	Acquisition system	33
8.2	BNC inputs CH1 to CH4	33
8.3	Trigger system	34
8.4	Physical	34
8.5	Power	34
8.6	Interface	34
8.7	I/O connectors	34
8.8	System requirements	35
8.9	Environmental conditions	35
8.10	Certifications and Compliances	35
8.11	Measure leads	35
8.12	Differential attenuators	36
8.13	Package contents	36

When working with electricity, no instrument can guarantee complete safety. It is the responsibility of the person who works with the instrument to operate it in a safe way. Maximum security is achieved by selecting the proper instruments and following safe working procedures. Safe working tips are given below:

- Always work according (local) regulations.
- Work on installations with voltages higher than 25 V_{AC} or 60 V_{DC} should only be performed by qualified personnel.
- Avoid working alone.
- Observe all indications on the ATS5004D before connecting any wiring
- Check the probes/test leads for damages. Do **not** use them if they are damaged
- Take care when measuring at voltages higher than 25 V_{AC} or 60 V_{DC}.
- Do not operate the equipment in an explosive atmosphere or in the presence of flammable gases or fumes.
- Do not use the equipment if it does not operate properly. Have the equipment inspected by qualified service personal. If necessary, return the equipment to GMTO for service and repair to ensure that safety features are maintained.



TiePie engineering
Koperslagersstraat 37
8601 WL Sneek
The Netherlands

EC Declaration of conformity

We declare, on our own responsibility, that the product

ATS5004D

for which this declaration is valid, is in compliance with

EN 55011:2009/A1:2010

IEC 61000-6-1/EN 61000-6-1:2007

EN 55022:2006/A1:2007

IEC 61000-6-3/EN 61000-6-3:2007

according the conditions of the EMC standard 2004/108/EC

and also with

Canada: ICES-001:2004

Australia/New Zealand: AS/NZS

Sneek, 1-11-2010

ir. A.P.W.M. Poelsma

Environmental considerations

This section provides information about the environmental impact of the ATS5004D.

ATS5004D end-of-life handling

Production of the ATS5004D required the extraction and use of natural resources. The equipment may contain substances that could be harmful to the environment or human health if improperly handled at the ATS5004D's end of life.

In order to avoid release of such substances into the environment and to reduce the use of natural resources, recycle the ATS5004D in an appropriate system that will ensure that most of the materials are reused or recycled appropriately.

The symbol shown below indicates that the ATS5004D complies with the European Union's requirements according to Directive 2002/96/EC on waste electrical and electronic equipment (WEEE).



Restriction of Hazardous Substances

The ATS5004D has been classified as Monitoring and Control equipment, and is outside the scope of the 2002/95/EC RoHS Directive.



Before using the ATS5004D first read chapter 1 about safety.

Many technicians investigate electrical signals. Though the measurement may not be electrical, the physical variable is often converted to an electrical signal, with a special transducer. Common transducers are accelerometers, pressure probes, current clamps and temperature probes. The advantages of converting the physical parameters to electrical signals are large, since many instruments for examining electrical signals are available.

The ATS5004D is a portable four channel measuring instrument with differential inputs. The native resolution is 12 bits, but user selectable resolutions of 14 and 16 bits are available too, with adjusted maximum sampling frequency:

resolution	Maximum sampling frequency
12 bit	50 MS/s
14 bit	3.125 MS/s
16 bit	195 kS/s

Table 3.1: Maximum sampling frequencies

With the accompanying software the ATS5004D can be used as an oscilloscope, a spectrum analyzer, a true RMS voltmeter or a transient recorder. All instruments measure by sampling the input signals, digitizing the values, process them, save them and display them.

3.1 Differential input

Most oscilloscopes are equipped with standard, single ended inputs, which are referenced to ground. This means that one side of the input is always connected to ground and the other side to the point of interest in the circuit under test.

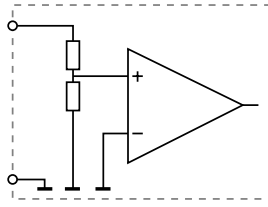


Figure 3.1: Single ended input

Therefore the voltage that is measured with an oscilloscope with standard, single ended inputs is always measured between that specific point and ground.

When the voltage is not referenced to ground, connecting a standard single ended oscilloscope input to the two points would create a short circuit between one of the points and ground, possibly damaging the circuit and the oscilloscope.

A safe way would be to measure the voltage at one of the two points, in reference to ground and at the other point, in reference to ground and then calculate the voltage difference between the two points. On most oscilloscopes this can be done by connecting one of the channels to one point and another channel to the other point and then use the math function $CH1 - CH2$ in the oscilloscope to display the actual voltage difference.

There are some disadvantages to this method:

- a short circuit to ground can be created when an input is wrongly connected
- to measure one signal, two channels are occupied
- by using two channels, the measurement error is increased, the errors made on each channel will be combined, resulting in a larger total measurement error
- The Common Mode Rejection Ratio (CMRR) of this method is relatively low. If both points have a relative high voltage, but the voltage difference between the two points is small, the voltage difference can only be measured in a high input range, resulting in a low resolution

A much better way is to use an oscilloscope with a differential input.

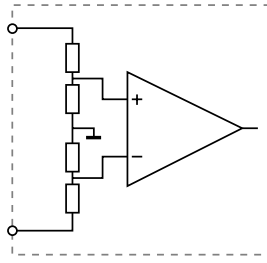


Figure 3.2: Differential input

A differential input is not referenced to ground, but both sides of the input are "floating". It is therefore possible to connect one side of the input to one point in the circuit and the other side of the input to the other point in the circuit and measure the voltage difference directly.

Advantages of a differential input:

- No risk of creating a short circuit to ground
- Only one channel is required to measure the signal
- More accurate measurements, since only one channel introduces a measurement error
- The CMRR of a differential input is high. If both points have a relative high voltage, but the voltage difference between the two points is small, the voltage difference can be measured in a low input range, resulting in a high resolution

3.1.1 Differential attenuators

To increase the input range of the ATS5004D, it comes with a differential 1:10 attenuator for each channel. This differential attenuator is specially designed to be used with the ATS5004D.



Figure 3.3: Differential attenuator

For a differential input, both sides of the input need to be attenuated.

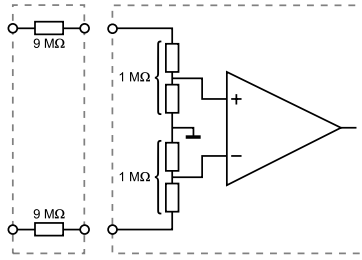


Figure 3.4: Differential input

Standard oscilloscope probes and attenuators only attenuate one side of the signal path. These are not suitable to be used with a differential input. Using these on a differential input will have a negative effect on the CMRR and will introduce measurement errors.

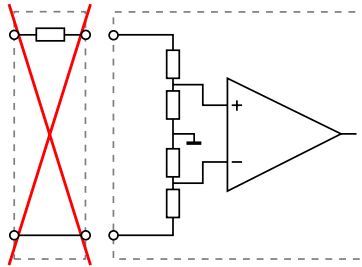


Figure 3.5: Differential input



The Differential Attenuator and the inputs of the Handy-scope HS4 Diff are differential, which means that the outside of the BNC's are not grounded, but carry life signals.

When using the attenuator, the following points have to be taken into consideration:

- do not connect other cables to the attenuator than the ones that are supplied with the instrument
- do not touch the metal parts of the BNC's when the attenuator is connected to the circuit under test, they can carry

a dangerous voltage. It will also influence the measurements and create measurement errors.

- do not connect the outside of the two BNC's of the attenuator to each other as this will short circuit a part of the internal circuit and will create measurement errors
- do not connect the outside of the BNC's of two or more attenuators that are connected to different channels of the Handyscope HS4 Diff to each other
- do not apply excessive mechanical force to the attenuator in any direction (e.g. pulling the cable, using the attenuator as handle to carry the ATS5004D, etc.)

3.1.2 Differential test lead

The ATS5004D comes with a special differential test lead. This test lead is specially designed to ensure a good CMRR.

The special heat resistant differential test lead provided with the ATS5004D is designed to be immune for noise from the surrounding environment.

3.2 Sampling

When sampling the input signal, samples are taken at fixed intervals. At these intervals, the size of the input signal is converted to a number. The accuracy of this number depends on the resolution of the instrument. The higher the resolution, the smaller the voltage steps in which the input range of the instrument is divided. The acquired numbers can be used for various purposes, e.g. to create a graph.

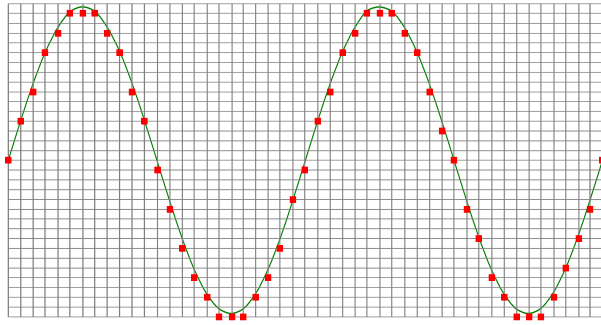


Figure 3.6: Sampling

The sine wave in figure 3.6 is sampled at the dot positions. By connecting the adjacent samples, the original signal can be reconstructed from the samples. You can see the result in figure 3.7.

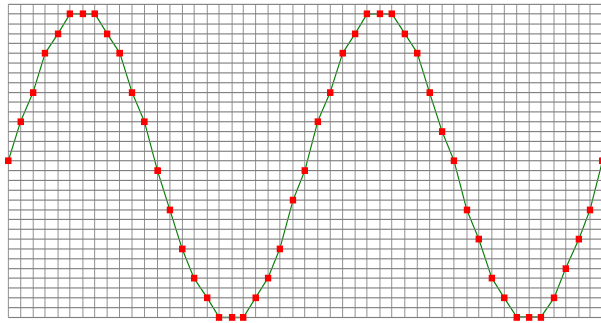


Figure 3.7: "connecting" the samples

3.3 Sample frequency

The rate at which the samples are taken is called the **sampling frequency**, the number of samples per second. A higher sampling frequency corresponds to a shorter interval between the samples. As is visible in figure 3.8, with a higher sampling frequency, the original signal can be reconstructed much better from the measured samples.

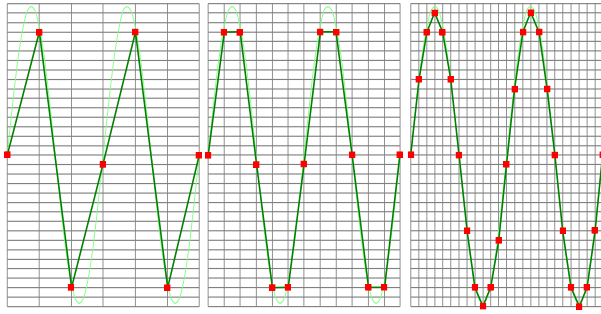


Figure 3.8: The effect of the sampling frequency

The sampling frequency must be higher than 2 times the highest frequency in the input signal. This is called the **Nyquist frequency**. Theoretically it is possible to reconstruct the input signal with more than 2 samples per period. In practice, 10 to 20 samples per period are recommended to be able to examine the signal thoroughly.

3.3.1 Aliasing

When sampling an analog signal with a certain sampling frequency, signals appear in the output with frequencies equal to the sum and difference of the signal frequency and multiples of the sampling frequency. For example, when the sampling frequency is 1000 Hz and the signal frequency is 1250 Hz, the following signal frequencies will be present in the output data:

Multiple of sampling frequency	1250 Hz signal	-1250 Hz signal
...		
-1000	$-1000 + 1250 = \mathbf{250}$	$-1000 - 1250 = -2250$
0	$0 + 1250 = 1250$	$0 - 1250 = -1250$
1000	$1000 + 1250 = 2250$	$1000 - 1250 = -250$
2000	$2000 + 1250 = 3250$	$2000 - 1250 = 750$
...		

Table 3.2: Aliasing

As stated before, when sampling a signal, only frequencies lower than half the sampling frequency can be reconstructed. In this case the sampling frequency is 1000 Hz, so we can only observe

signals with a frequency ranging from 0 to 500 Hz. This means that from the resulting frequencies in the table, we can only see the 250 Hz signal in the sampled data. This signal is called an **alias** of the original signal.

If the sampling frequency is lower than twice the frequency of the input signal, **aliasing** will occur. The following illustration shows what happens.

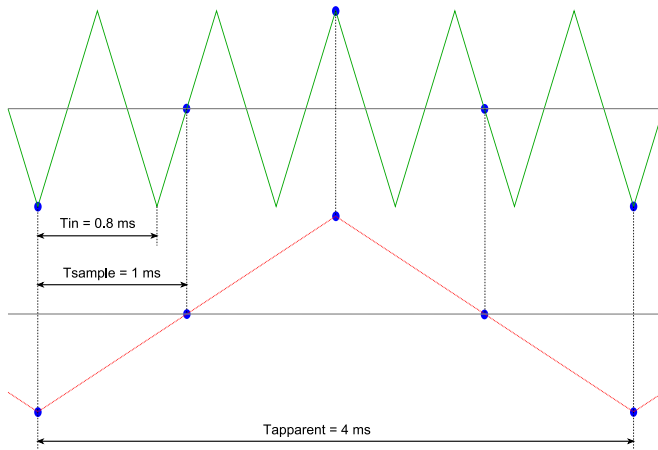


Figure 3.9: Aliasing

In figure 3.9, the green input signal (top) is a triangular signal with a frequency of 1.25 kHz. The signal is sampled with a frequency of 1 kHz. The corresponding sampling interval is $1/1000\text{Hz} = 1\text{ms}$. The positions at which the signal is sampled are depicted with the blue dots. The red dotted signal (bottom) is the result of the reconstruction. The period time of this triangular signal appears to be 4 ms, which corresponds to an apparent frequency (alias) of 250 Hz (1.25 kHz - 1 kHz).



To avoid aliasing, always start measuring at the highest sampling frequency and lower the sampling frequency if required.

3.4 Digitizing

When digitizing the samples, the voltage at each sample time is converted to a number. This is done by comparing the voltage with a number of levels. The resulting number is the number corresponding to the level that is closest to the voltage. The number of levels is determined by the resolution, according to the following relation: $LevelCount = 2^{Resolution}$.

The higher the **resolution**, the more levels are available and the more accurate the input signal can be reconstructed. In figure 3.10, the same signal is digitized, using two different amounts of levels: 16 (4-bit) and 64 (6-bit).

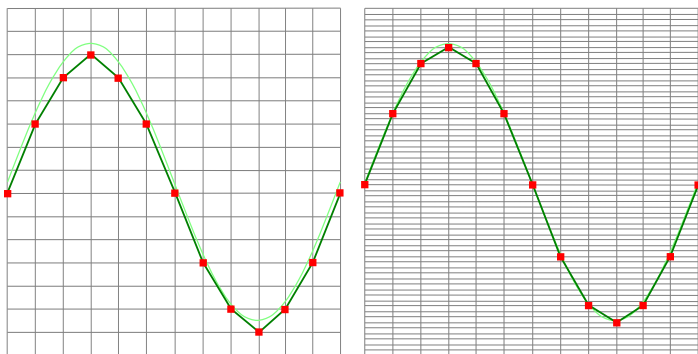


Figure 3.10: The effect of the resolution

The ATS5004D measures at e.g. 12 bit resolution ($2^{12}=4096$ levels). The smallest detectable voltage step depends on the input range. This voltage can be calculated as:

$$VoltageStep = FullInputRange / LevelCount$$

For example, the 200 mV range ranges from -200 mV to +200 mV, therefore the full range is 400 mV. This results in a smallest detectable voltage step of $0.400V / 4096 = 97.65 \mu V$.

3.5 Signal coupling

The ATS5004D has two different settings for the signal coupling: AC and DC. In the setting DC, the signal is directly coupled to the

input circuit. All signal components available in the input signal will arrive at the input circuit and will be measured.

In the setting AC, a capacitor will be placed between the input connector and the input circuit. This capacitor will block all DC components of the input signal and let all AC components pass through. This can be used to remove a large DC component of the input signal, to be able to measure a small AC component at high resolution.



When measuring DC signals, make sure to set the signal coupling of the input to DC.



Before connecting the ATS5004D to the computer, the drivers need to be installed.

4.1 Introduction

To operate a ATS5004D, a driver is required to interface between the measurement software and the instrument. This driver takes care of the low level communication between the computer and the instrument, through USB. When the driver is not installed, or an old, no longer compatible version of the driver is installed, the software will not be able to operate the ATS5004D properly or even detect it at all.

The installation of the USB driver is done in a few steps. Firstly, the driver has to be pre-installed by the driver setup program. This makes sure that all required files are located where Windows can find them. When the instrument is plugged in, Windows will detect new hardware and install the required drivers.

4.2 Where to find the driver setup

The driver setup program and measurement software can be found in the download section on TiePie engineering's website and on the CD-ROM that came with the instrument. It is recommended to install the latest version of the software and USB driver from the website. This will guarantee the latest features are included.

4.3 Executing the installation utility

To start the driver installation, execute the downloaded driver setup program, or the one on the CD-ROM that came with the instrument. The driver install utility can be used for a first time

installation of a driver on a system and also to update an existing driver.

The screen shots in this description may differ from the ones displayed on your computer, depending on the Windows version.

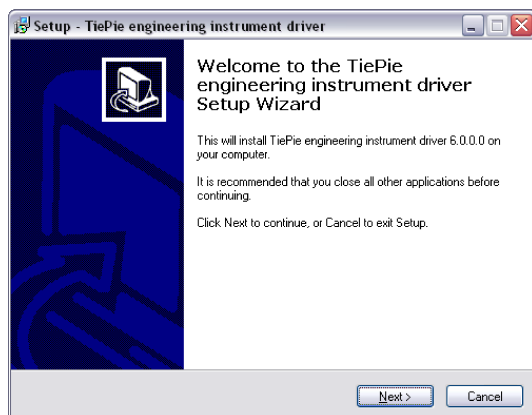


Figure 4.1: Driver install: step 1

When drivers were already installed, the install utility will remove them before installing the new driver. To remove the old driver successfully, **it is essential** that the ATS5004D is disconnected from the computer prior to starting the driver install utility. When the ATS5004D is used with an external power supply, this must be disconnected too.

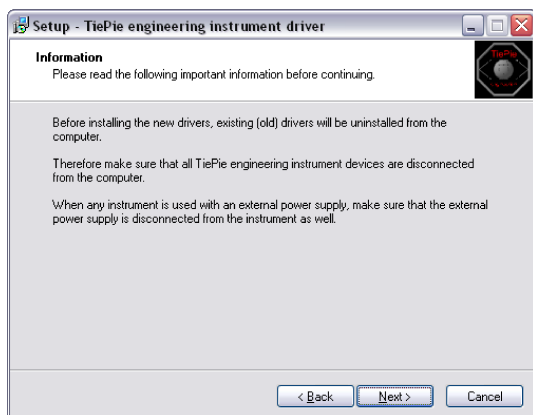


Figure 4.2: Driver install: step 2

When the instrument is still connected, the driver install utility will recognize it and report this. You will be asked to continue anyway.

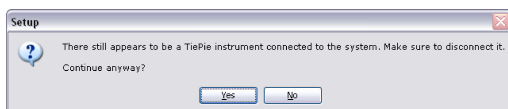


Figure 4.3: Driver install: Instrument is still connected

Clicking "**No**" will bring back the previous screen. The instrument should now be disconnected. Then the removal of the existing driver can be continued by clicking "**Next**".

Clicking "**Yes**" will ignore the fact that the instrument is still connected and continue removal of the old driver. This option is **not** recommended, as removal may fail, after which installation of the new driver may fail as well.

When no existing driver was found or the existing driver is removed, the location for the pre-installation of the new driver can be selected.

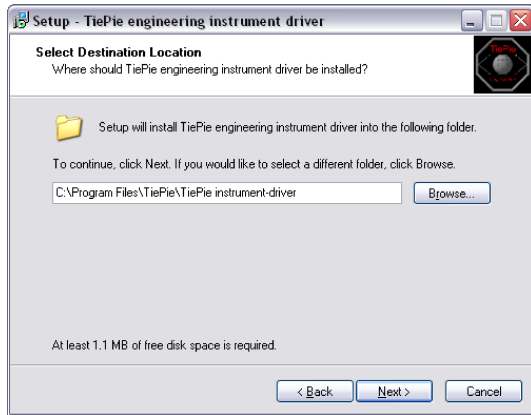


Figure 4.4: Driver install: step 3

On Windows XP and newer, the installation may inform about the drivers not being "Windows Logo Tested". The driver is not causing any danger for your system and can be safely installed. Please ignore this warning and continue the installation.

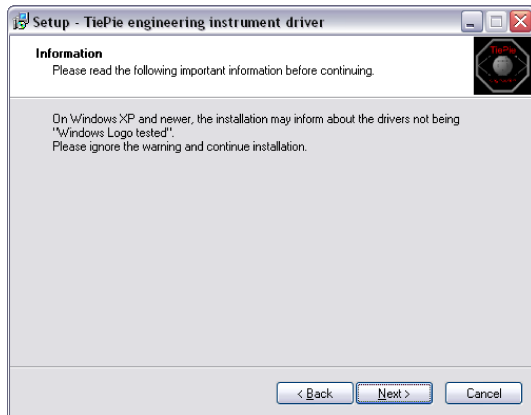


Figure 4.5: Driver install: step 4

The driver install utility now has enough information and can install the drivers. Clicking "**Install**" will remove existing drivers and install the new driver. A *remove entry* for the new driver is added to the software applet in the Windows control panel.

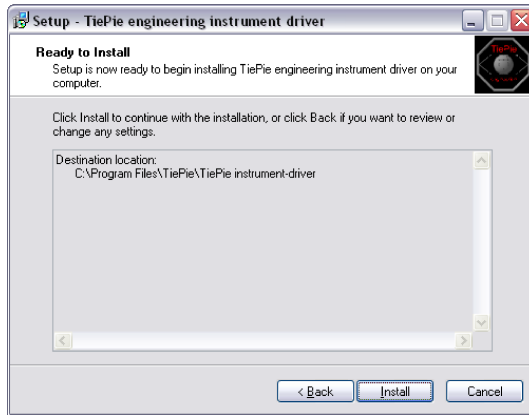


Figure 4.6: Driver install: step 5

As mentioned, Windows XP SP2 and newer may warn for the USB drivers not being Windows Logo tested. Please ignore this warning and continue anyway.



Figure 4.7: Driver install: Ignore warning and continue

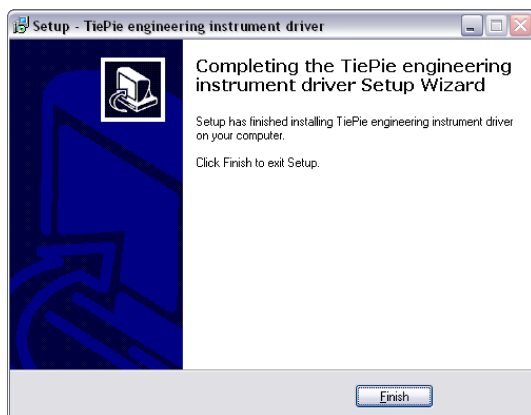


Figure 4.8: Driver install: Finished



Drivers have to be installed before the ATS5004D is connected to the computer for the first time. See chapter 4 for more information.

5.1 Power the instrument

The ATS5004D is powered by the USB, no external power supply is required. Only connect the ATS5004D to a bus powered USB port, otherwise it may not get enough power to operate properly.

5.1.1 External power

In certain cases, it can be that the ATS5004D cannot get enough power from the USB port.

When a ATS5004D is connected to a USB port, the hardware will be powered, resulting in an inrush current, which is higher than the nominal current. After the inrush current, the current will stabilize at the nominal current.

USB ports have a maximum limit for both the inrush current peak and the nominal current. When either of them is exceeded, the USB port will be switched off. As a result, the connection to the ATS5004D will be lost.

Most USB ports can supply enough current for the ATS5004D to work without an external power supply, but this is not always the case. Some (battery operated) portable computers or (bus powered) USB hubs do not supply enough current. The exact value at which the power is switched off, varies per USB controller, so it is possible that the ATS5004D functions properly on one computer, but does not on another.

In order to power the ATS5004D externally, an external power input is provided for. It is located at the rear of the ATS5004D. Refer to paragraph 7.1 for specifications of the external power input.

5.2 Connect the instrument to the computer

After the new driver has been pre-installed (see chapter 4), the ATS5004D can be connected to the computer. When the ATS5004D is connected to a USB port of the computer, Windows will report new hardware. The Found New Hardware Wizard will appear.

Depending on the Windows version, the New Hardware Wizard will show a number of screens in which it will ask for information regarding the drivers of the newly found hardware. The appearance of the dialogs will differ for each Windows version and might be different on the computer where the ATS5004D is installed.



The driver consists of two parts which are installed separately.

Once the first part is installed, the installation of the second part will start automatically. Installation of the second part is identical to the first part, therefore they are not described individually here.

5.2.1 Found New Hardware Wizard

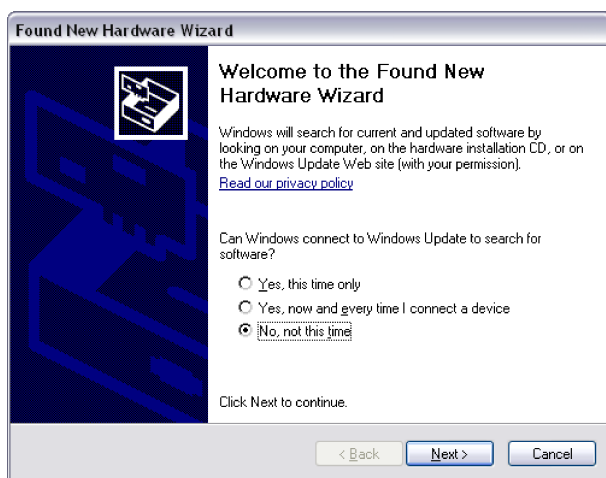


Figure 5.1: Hardware install: step 1

This window will only be shown in Windows XP SP2 or newer. No drivers for the ATS5004D can be found on the Windows Update Web site, so select *"No, not this time"* and click *"Next"*.

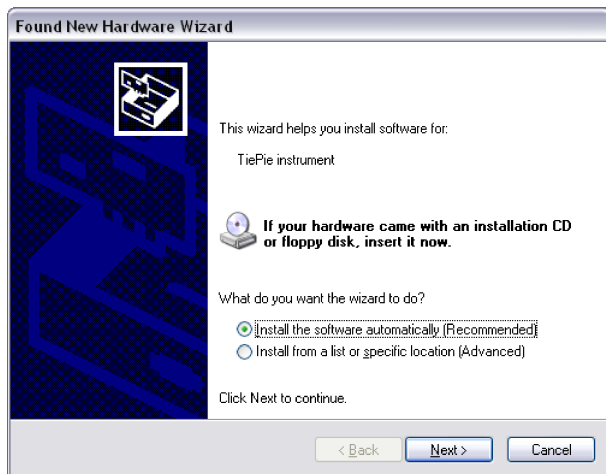


Figure 5.2: Hardware install: step 2

Since the drivers are already pre-installed on the computer, Windows will be able to find them automatically. Select *"Install the software automatically"* and click *"Next"*.

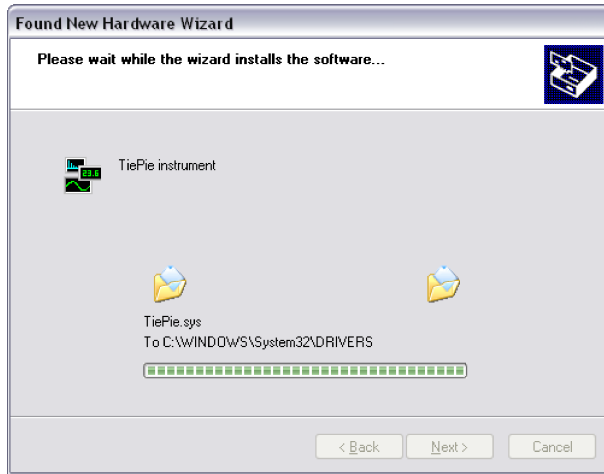


Figure 5.3: Hardware install: step 3

The New Hardware wizard will now copy the required files to their destination.

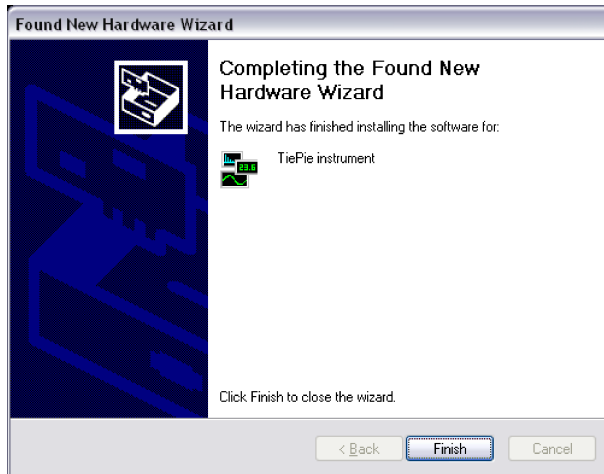


Figure 5.4: Hardware install: step 4

The first part of the new driver is now installed. Click "**Finish**" to close the wizard and start installation of the second part, which follows identical steps.

Once the second part of the driver is installed, measurement software can be installed and the ATS5004D can be used.

5.3 Plug into a different USB port

When the ATS5004D is plugged into a different USB port, some Windows versions will treat the ATS5004D as different hardware and will ask to install the drivers again. This is controlled by Microsoft Windows and is not caused by GMT0.

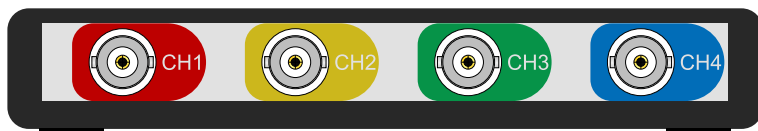


Figure 6.1: Front panel

6.1 Channel input connectors

The CH1 – CH4 BNC connectors are the main inputs of the acquisition system. The isolated BNC connectors are not connected to the ground of the ATS5004D.

6.2 Power indicator

A power indicator is situated at the top cover of the instrument. It is lit when the ATS5004D is powered.

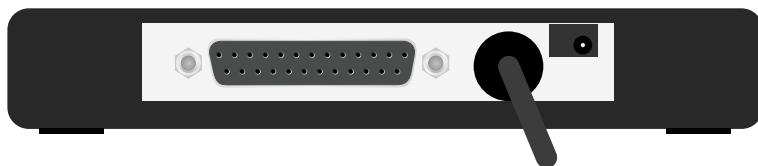


Figure 7.1: Rear panel

7.1 Power

The ATS5004D is powered through the USB. If the USB cannot deliver enough power, it is possible to power the instrument externally. The ATS5004D has two external power inputs located at the rear of the instrument: the dedicated power input and a pin of the extension connector.

The specifications of the dedicated power connector are:



Pin	Dimension	Description
Center pin	Ø1.3 mm	ground
Outside bushing	Ø3.5 mm	positive

Figure 7.2: Power connector

Besides the external power input, it is also possible to power the instrument through the extension connector, the 25 pin D-sub connector at the rear of the instrument. The power has to be applied to pin 3 of the extension connector. Pin 4 can be used as ground.

The following minimum and maximum voltages apply to both power inputs:

Minimum	Maximum
4.5 V _{DC}	14 V _{DC}

Table 7.1: Maximum voltages

Note that the externally applied voltage should be higher than the USB voltage to relieve the USB port.

7.1.1 USB power cable

The ATS5004D is delivered with a special USB external power cable.



Figure 7.3: USB power cable

One end of this cable can be connected to a second USB port on the computer, the other end can be plugged in the external power input at the rear of the instrument. The power for the instrument will be taken from two USB ports of the computer.



The outside of the external power connector is connected to +5 V. In order to avoid shortage, first connect the cable to the ATS5004D and then to the USB port.

7.1.2 Power adapter

In case a second USB port is not available, or the computer still can't provide enough power for the instrument, an external power adapter can be used. When using an external power adapter, make sure that:

- the polarity is set correctly
- the voltage is set to a valid value for the instrument and higher than the USB voltage
- the adapter can supply enough current (preferably >1 A)
- the plug has the correct dimensions for the external power input of the instrument

7.2 USB

The ATS5004D is equipped with a USB 2.0 High speed (480 Mbit/s) interface with a fixed cable with type A plug. It will also work on a computer with a USB 1.1 interface, but will then operate at 12 Mbit/s.

7.3 Extension Connector

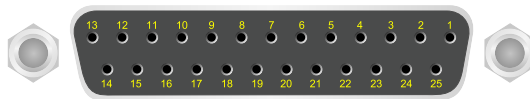


Figure 7.4: Extension connector

To connect to the ATS5004D a 25 pin female D-sub connector is available, containing the following signals:

Pin	Description	Pin	Description
1	Ground	14	Ground
2	Reserved	15	Ground
3	External Power in DC	16	Reserved
4	Ground	17	Ground
5	+5V out, 10 mA max.	18	Reserved
6	Ext. sampling clock in (TTL)	19	Reserved
7	Ground	20	Reserved
8	Ext. trigger in (TTL)	21	Reserved
9	Data OK out (TTL)	22	Ground
10	Ground	23	I ² C SDA
11	Trigger out (TTL)	24	I ² C SCL
12	Reserved	25	Ground
13	Ext. sampling clock out (TTL)		

Table 7.2: Pin description Extension connector

All TTL signals are 3.3 V TTL signals which are 5 V tolerant, so they can be connected to 5 V TTL systems.

Pins 9, 11, 12, 13 are open collector outputs. Connect a pull-up resistor of 1 kOhm to pin 5 when using one of these signals.

8.1 Acquisition system

Number of input channels	4 analog
CH1, CH2, CH3, CH4	isolated BNC
Maximum sampling rate	depending on resolution
12 bit	50 MS/s
14 bit	3.125 MS/s
16 bit	195.3 kS/s
Maximum streaming rate	depending on resolution
12 bit	500 kS/s
14 bit	480.8 kS/s
16 bit	195.3 kS/s
Sampling source	internal quartz, external
Internal	Quartz
Accuracy	$\pm 0.01\%$
Stability	± 100 ppm over -40°C to $+85^{\circ}\text{C}$
External	On extension connector
Voltage	3.3 V TTL, 5 V TTL tolerant
Frequency range	95 MHz to 105 MHz
Memory	128 kSamples per channel

8.2 BNC inputs CH1 to CH4

Type	Differential
Resolution	12, 14, 16 bit user selectable
Accuracy	0.3% of full scale ± 1 LSB
Range	200 mV to 80 V full scale
Coupling	AC/DC
Impedance	1 M Ω / 30 pF
Maximum voltage	200 V (DC + AC peak <10 kHz)
Maximum voltage with 1:10 attenuator	300 V (DC + AC peak <10 kHz)
Maximum common mode voltage	200 mV to 800 mV ranges : 2 V 2 V to 8 V ranges : 20 V 20 V to 80 V ranges : 200 V
Common Mode Rejection Ratio	-48 dB
Channel Isolation	500 V
Channel Separation	-80 dB
Bandwidth (-3dB)	50 MHz
AC coupling cut off frequency (-3dB)	± 1.5 Hz

8.3 Trigger system

System	digital, 2 levels
Source	CH1, CH2, CH3, CH4, digital external, AND, OR
Trigger modes	rising slope, falling slope, inside window, outside window
Level adjustment	0 to 100% of full scale
Hysteresis adjustment	0 to 100% of full scale
Resolution	0.024 % (12 bits)
Pre trigger	0 to 128 ksamples (0 to 100%, one sample resolution)
Post trigger	0 to 128 ksamples (0 to 100%, one sample resolution)
Trigger hold-off	0 to 1 MSamples, 1 sample resolution
Digital external trigger	
Input	extension connector
Range	0 to 5 V (TTL)
Coupling	DC

8.4 Physical

Instrument height	25 mm / 1.0"
Instrument length	170 mm / 6.7"
Instrument width	140 mm / 5.2"
Weight	480 gram / 17 ounce
USB cord length	1.8 m / 70"

8.5 Power

Input	from USB or external input
Consumption	500 mA max

8.6 Interface

Interface	USB 2.0 High Speed (480 Mbit/s) (USB 1.1 Full Speed (12 Mbit/s) and USB 3.0 compatible)
-----------	---

8.7 I/O connectors

CH1 .. CH4	BNC
Power	3.5 mm power socket
Extension connector	D-sub 25 pins female
USB	Fixed cable with type A plug

8.8 System requirements

PC I/O connection	USB 2.0 High Speed (480 Mbit/s) (USB 1.1 and USB 3.0 compatible)
Operating System	Windows XP / Vista / 7 / 8 / 10

8.9 Environmental conditions

Operating	
Ambient temperature	0°C to 55°C
Relative humidity	10 to 90% non condensing
Storage	
Ambient temperature	-20°C to 70°C
Relative humidity	5 to 95% non condensing

8.10 Certifications and Compliances

CE mark compliance	Yes
RoHS	Yes

8.11 Measure leads

Model	TP-C812B
Type	differential
Connectors	
Instrument side	isolated female BNC connector
Test point side	red and black 4 mm shrouded banana jacks
Bandwidth	7 MHz
Safety	CAT III, 1000 V, double isolated
Dimensions	
Total length	2000 mm
Length to split	800 mm
Length individual ends	1200 mm
Weight	80 g
Color	black
Heat resistant	yes
Certification and compliances	
CE conformity	yes
RoHS	yes
Accessories	
Color coding rings	5 x 3 rings, various colors
Suitable instrument	ATS5004D

8.12 Differential attenuators

Model	TP-DA10
Attenuation	10 x, differential
Bandwidth	25 MHz
Input impedance	10 M Ω // 15 pF
Maximum input voltage	300 V
Input connector	Female BNC
Output connector	Male BNC
Dimensions	
Length	79 mm
Diameter	19 mm
Weight	30 g
Suitable instrument	ATS5004D

8.13 Package contents

Instrument	ATS5004D
Test lead	4 x TP-C812B low noise differential test leads with 4 mm banana jacks
Accessories	4 x TP-DA10 1:10 differential attenuators USB power cable
Software	Windows XP / Vista / 7 / 8 / 10
Drivers	Windows XP / Vista / 7 / 8 / 10
Manual	Instrument manual and software user's manual

